

Technical Specifications

Output	Up to 22 J/Cm ²
Treatment Filters	HR: 530-1200 nm SR: 570-1200 nm AC: 400-1200 nm
Pulse Duration	25 msec
Pulse Repetition Rate	0.33Hz (every 2-3 seconds)
Spot Size	15 X 50 mm (7.5cm ²)
Dimensions	45cm (L) x 38cm (W) x 27cm (H)
Weight	17 kg / 37lbs
Electrical	220V/50Hz or 110V/60Hz - 800 V/A

About Viora

Viora is a market leader in the development of technologically advanced medical aesthetic systems for physicians and aesthetic clinicians worldwide. Dedicated to the research and development of safe and effective aesthetic devices, Viora provides the highest standard of professional grade solutions while keeping your business in mind.

Viora offers a wide range of applications including body contouring, cellulite reduction, and skin tightening procedures as well as phototherapy and anti-aging treatments. All of the products are designed by a team of veteran medical device developers in close cooperation with the company's advisory board including medical aesthetics practitioners, plastic surgeons, and dermatologists. Viora's aesthetic devices are clinically tested and endorsed by physicians and aesthetic specialists worldwide.

The Trios™ system is manufactured under ISO 13485.

Manufacturer's Address

Viora Ltd.
3 Maskit Street
46733
Herzliya, Israel
Tel: +972 9955 1344
Fax: +972 9955 1345
Website: www.vioramed.com
E-mail: info@vioramed.com

Manufacturer's Authorized Representative

CEPartner4U BV
Esdoornlaan 13
3951 DB Maarn
The Netherlands
Tel: +31 (0) 343 442 524
Fax: +31 (0) 343 442 162

FDA cleared for marketing

This device complies with the European Medical Device Directive 93/42/EEC



TRIOS

Smooth, Radiant Skin in a Flash



A d v a n c e d P h o t o t h e r a p y i n a N e w L i g h t



It All Begins Here

Created in the heart of the world's IPL industry, Trios™ was developed by a research and development team with over two decades of engineering experience in IPL based systems. Trios™ is a phototherapy system that delivers the most sought after IPL applications in the aesthetic beauty market from long-term hair removal to acne clearance.

As an award-winning company, Viora takes pride in delivering high-end, premium quality systems that affirm our commitment and partnership to promoting your business needs and services.

The Ability to Offer More

Trios™ empowers your business by providing you with the industry's number one phototherapy procedure, hair removal. And now, you can offer even more in phototherapy with skin rejuvenation and acne clearance treatments.

Trios™ is a modular system with versatile applications that enable you to expand your treatment portfolio and enlarge your customer database by offering a suite of aesthetic solutions:


- Long-term hair removal for normal hair growth & Hirsutism
- Rejuvenates skin, reducing vascular and pigmented lesions
- Smooths fine lines and wrinkles
- Balances skin pigmentation caused by sun damage & age spots
- Acne clearance

Results, Right on Target


Based on a process called selective photothermolysis, Trios™ emits pulsed light energy in one single flash. The light emitted is absorbed by targeted chromophores (light sensitive cells) and encourages a specific biological process in order to achieve the desired clinical result.




Hair Removal

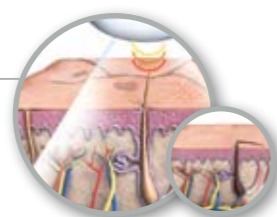

Skin Rejuvenation

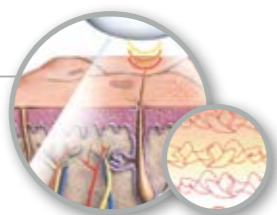

Acne Clearance

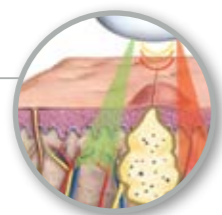

A Host of Advantages

Trios™ was designed to ensure practitioners are equipped with a system that delivers uncompromising performance, safety and efficiency on all levels.

To this end, Trios™ offers advantages that are second to none on the market:

- Actual 22 J/cm² of power at skin level for effective treatments and maximum results.
- Highly accurate program settings enable you to tailor the best treatment for your patients.
- Highly reliable device; manufactured to the highest standards and requires minimal setup & maintenance.
- Safety in mind: advanced computerized management control to ensure system performance and risk-free treatments.
- Visible results after the first treatment.
- Fast, comfortable treatments.
- Digital, user-friendly interface.
- Unique, ergonomic hand piece design.
- Portable device allows users to move between rooms and salon locations.



The Profitable Formula

Large spot size area (7.5cm²)

- Reduced treatment time
- Increased accuracy and efficacy
- Less pulses required per treatment

Long lamp life (100,000 pulses)

- Cost efficient

Fast repetition rate

- Faster treatment times
- Perform more treatment procedures per day

Multiple treatment applications

- Attract new customers
- Increase treatment offering to existing customers

